

# ClassicStar® NIV Full-Face Mask Consumables and Accessories


ClassicStar NIV full-face mask is designed for non-invasive ventilation therapy. Because the mask covers the nose and mouth, it provides effective therapy, even in case of mouth breathing. This portfolio supports premium comfort with anatomical fit and real seal.

#### ClassicStar® SE (with Standard Elbow)

The blue elbow indicates that the mask has neither an Anti-Asphyxia Valve nor exhalation ports.

#### ClassicStar® AAV (with Anti-Asphyxia Valve)

The white transparent elbow includes an Anti-Asphyxia Valve that ensures convenient access to room air in the event of reduced pressure (< 3 cm  $\rm H_2O$ ). Exhalation ports support  $\rm CO_2$  elimination during non-invasive ventilation.

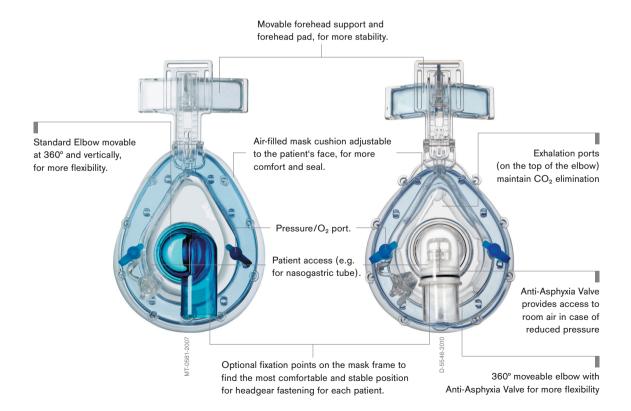

# **Benefits**

### ClassicStar® NIV Full-Face Mask

The ClassicStar NIV full-face mask is for single patient use.

The Dräger Sizing Gauge helps you to choose the right mask size - there is no need to open the mask package.

- The ClassicStar NIV full-face mask has an air-filled mask cushion, which can be adjusted individually
  to the patient's face. The mask cushion can be inflated or deflated by means of a pump ball. Due to its
  adaptability, the ClassicStar mask provides increased comfort and effective seal.
- The ClassicStar mask cushion allocates pressure over a larger surface area, which helps to avoid skin irritation.
- The forehead support distributes pressure from the bridge of the nose. The forehead pad and forehead support are adjustable to match the patient's forehead curve.
- The 6-point headgear provides a stable and secure fixation, thus reducing mask displacement and pressure. Safety release is integrated into the headgear for quick removal of the mask.
- Optional fixation points on the mask frame allow for the choice of the most convenient and stable headgear fastening position for each patient.
- The ClassicStar swivel is moveable at 360°. In addition, the swivel can be moved vertically (SE version only), providing more flexibility.

# Ordering Information

ClassicStar SE	Packing unit	Order number
ClassicStar SE, full-face mask for	1 pc.	MP01573
non-invasive ventilation with Standard Elbow,		
size S		
ClassicStar SE, full-face mask for	1 pc.	MP01574
non-invasive ventilation with Standard Elbow,		
size M		
ClassicStar SE, full-face mask for	1 pc.	MP01575
non-invasive ventilation with Standard Elbow,		
size L		
ClassicStar AAV		
ClassicStar AAV, full-face mask for	1 pc.	MP01570
non-invasive ventilation with		
Anti-Asphyxia Valve, size S		
ClassicStar AAV, full-face mask for	1 pc.	MP01571
non-invasive ventilation with		
Anti-Asphyxia Valve, size M		
ClassicStar AAV, full-face mask	1 pc.	MP01572
for non-invasive ventilation with		
Anti-Asphyxia Valve, size L		
Two-sided pump ball for air-managment	1 pc.	MP01590
Sizing Gauge for ClassicStar SE and AAV	3 pcs.	MP01592
masks for non-invasive ventilation		

## Notes

### CORPORATE HEADQUARTERS

Drägerwerk AG & Co. KGaA Moislinger Allee 53–55 23558 Lübeck, Germany www.draeger.com

#### Manufacturer:

Dräger Medical GmbH Moislinger Allee 53-55 23558 Lübeck, Germany

## As of August 2015

Dräger Medical GmbH changes to Drägerwerk AG & Co. KGaA

Locate your Regional Sales Representative at: www.draeger.com/contact



## REGION EUROPE CENTRAL AND EUROPE NORTH

Dräger Medical GmbH Moislinger Allee 53-55 23558 Lübeck, Germany Tel +49 451 882 0 Fax +49 451 882 2080 info@draeger.com

#### REGION EUROPE SOUTH

Dräger Médical S.A.S.
Parc de Haute Technologie
d'Antony 2
25, rue Georges Besse
92182 Antony Cedex, France
Tel +33 1 46 11 56 00
Fax +33 1 40 96 97 20
dlmfr-contact@draeger.com

## REGION MIDDLE EAST, AFRICA

Dräger Medical GmbH Branch Office P.O. Box 505108 Dubai, United Arab Emirates Tel +971 4 4294 600 Fax +971 4 4294 699 contactuae@draeger.com

#### REGION ASIA / PACIFIC

Draeger Medical South East Asia Pte Ltd. 25 International Business Park #04-27/29 German Centre Singapore 609916, Singapore Tel +65 6572 4388 Fax +65 6572 4399 asia.pacific@draeger.com REGION CENTRAL AND SOUTH AMERICA

Dräger Panama Comercial S. de R.L. Complejo Business Park, V tower, 10th floor Panama City Tel +507 377 9100 Fax +507 377 9130 contactcsa@draeger.com

The quality management system at Dräger Medical GmbH is certified according to ISO 13485, ISO 9001 and Annex II.3 of Directive 93/42/EEC (Medical devices).